



GENOTYPING PROTOCOL

Investigator: Mark Magnuson

Genome Edit: Zfp329 KO

Allele name: *Zfp329^{em1Mgn}*

Primers:

Zfp329Fwd: ATGGAGGGATTACAAAGAGAGG

Zfp329Rev: GGCAACTATGTAAGGTTTGGTC

Predicted PCR Product: WT = 271 bp, KO = 263 bp

Validated PCR protocol:

20 μ L PCR reaction

5x Phusion HF buffer = 4.0 μ L

10 mM dNTPs = 0.4 μ L

10 μ M Zfp329Fwd: = 1.0 μ L

10 μ M Zfp329Rev: = 1.0 μ L

Phusion Polymerase = 0.2 μ L

Nuclease-free water = 12.4 μ L

Genomic DNA (about 50 ng) = 1.0 μ L

PCR program

1. 98°C, 2 min

2. 98°C, 10 sec

3. 64°C, 15 sec


4. 72°C, 20 sec

5. Go to 2, 35X


6. 72°C, 7 min

7. 4°C, ∞

Run on agarose gel $\geq 2.5\%$ (improved resolution with longer runs and higher agarose gel percentage)

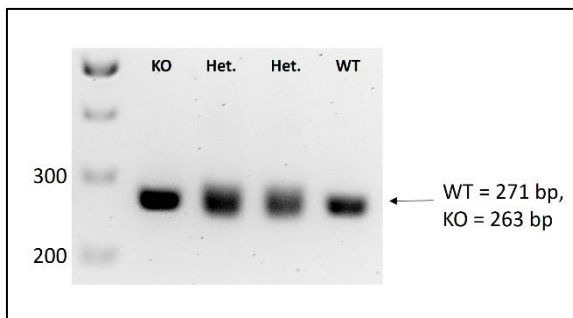

Note that differentiating between the WT and KO is challenging on an agarose gel for this assay. Higher percentage agarose gels may allow for differentiation, but untested.

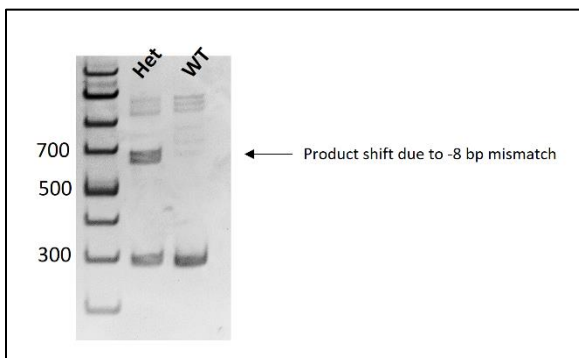

It is recommended to either sequence validate the samples or spike in a known WT PCR product at a 50:50 ratio in all the non-heterozygous samples, heat to 95°C for 10 minutes, then allow to cool to room temperature. Run the randomly annealed PCR product into a 10% TBE gel and identify KO samples by presence of a heterozygous doublet.