




Vanderbilt Genome Editing Resource

GENOTYPING PROTOCOL: *Sox17Rr* -391 bp

Investigator: Mark Magnuson
Genome Edit: Sox17Rr -391 bp
Allele Name: *Rr*_{em1Mgn} (Rr ID number is assigned upon publication)

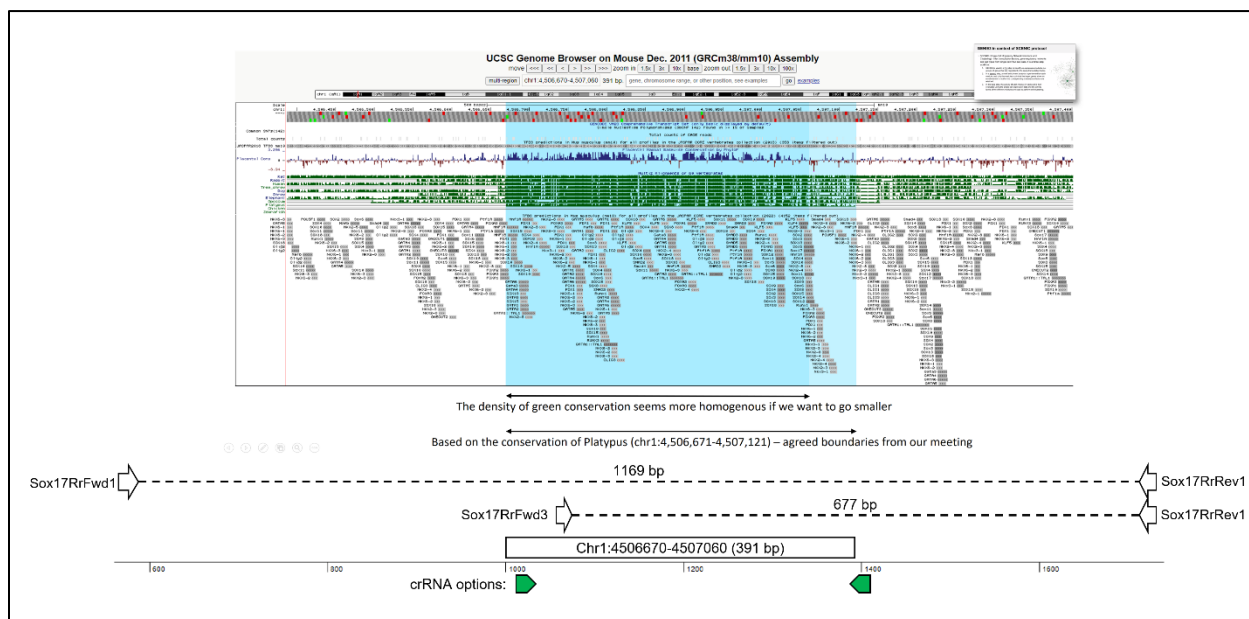

Primers:


Primer name	Sequence (5'-3')
Sox17RrFwd1	CCTGGTTATGAACTCTGAGAACAC
Sox17RrFwd3	AGGTTGCTGCTTGTAAGAGC
Sox17RrRev1	GAGAAGACCTCTGTGTACACC

PCR product sizes:

Allele	Primer set Sox17RrFwd1 + Rev1	Primer set Sox17RrFwd3 + Rev1
WT	1169 bp	677 bp
Sox17Rr – 391 bp	784 bp	none

PCR Protocol:

Polymerase	Phusion High-Fidelity DNA Polymerase NEB, catalog #M0530S
Cycling Conditions	Initial denaturation, 98°C, 30 seconds = 1 cycle Denaturation, 98°C, 10 seconds = 35 cycles Annealing, 62°C, 10 seconds Extension, 72°C, 25-35 seconds Final extension, 72°C, 2 minutes = 1 cycle Hold, 4°C, indefinitely
Final reaction concentrations (NEB protocol recommendations)	1x Phusion HF buffer 200 uM dNTPs 0.5 uM primer 1 0.5 uM primer 2 1 unit per 50 ul PCR Phusion DNA polymerase Variable template DNA input Nuclease-free water to volume

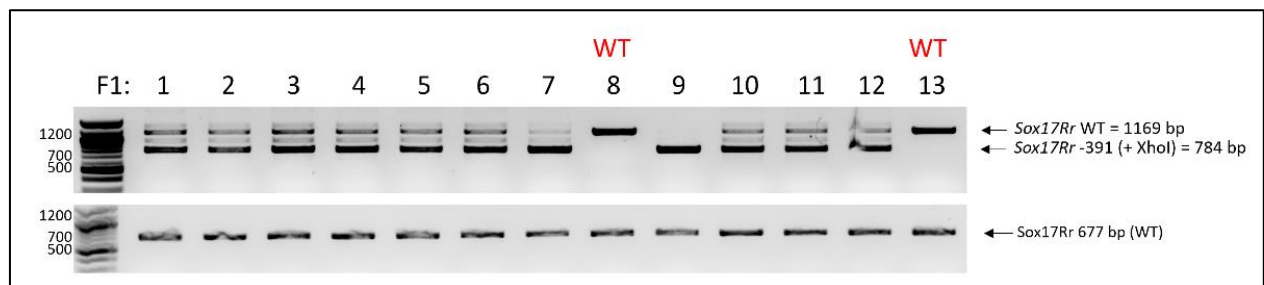Example image:**Commercial genotyping option:**

A TransnetYX qPCR assay is in progress for your strain.